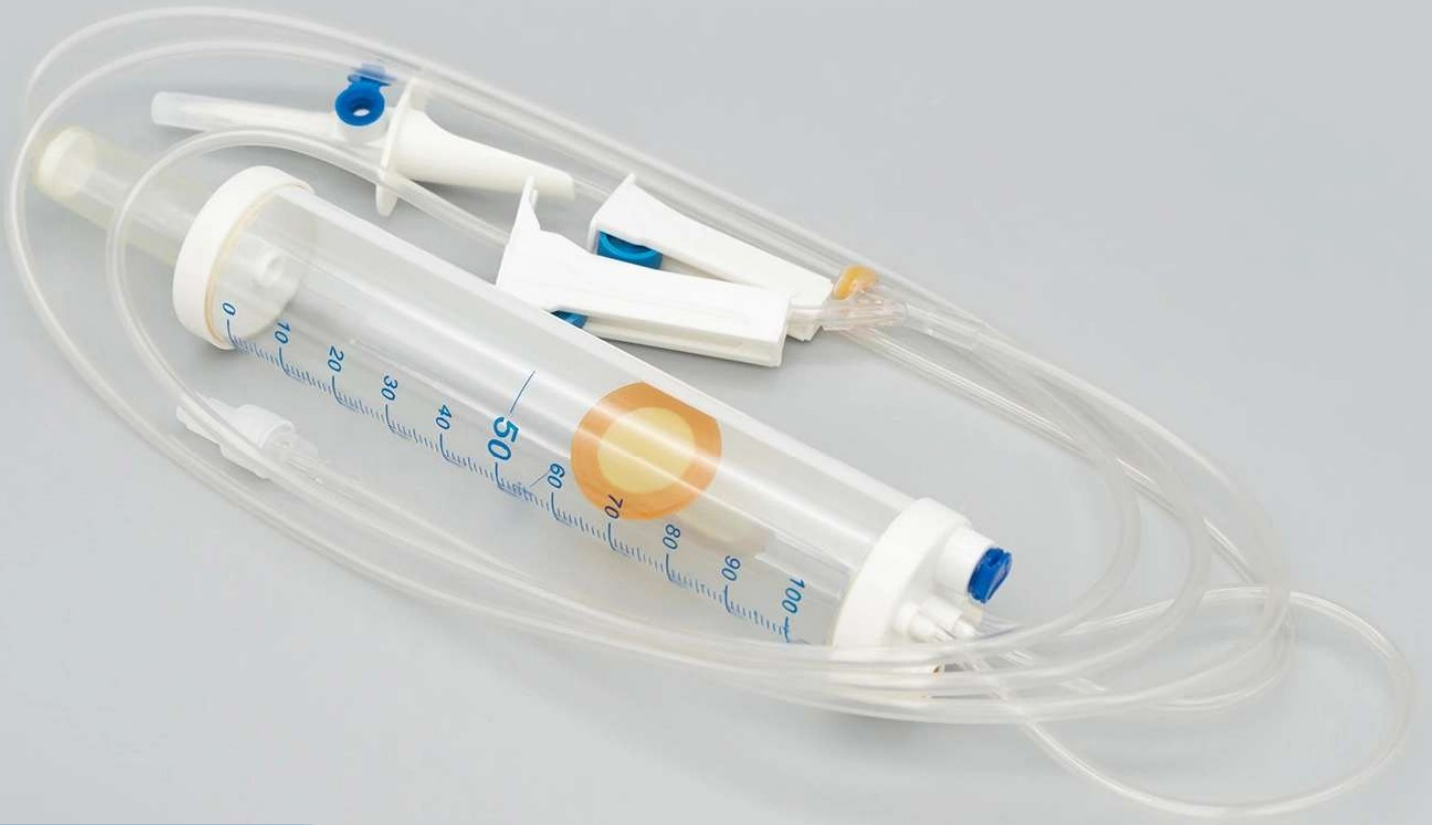




Infusion Burette Set

 Building a Healthier
Tomorrow

RAMY

1st Industrial City | P.O. Box 7917 |
Road 22 | Cross 11 | Dammam-31472 |
Kingdom of Saudi Arabia

Al Shifa Medical Co. located at Eastern Region 1st Industrial City, Dammam Kingdom of Saudi Arabia since **1980's**. Al Shifa supplies with a vast range of Medical Devices in the Middle East as well as exports to the different countries of Middle East, Europe, Asia, and Africa.

Al-Shifa supplies different types and sizes of Infusion Burette Set commonly used for patient who needs IV Therapy administration.

The design of Infusion Burette Set provides a vertical transparent cylinder shape with a volumetric graduation marking on its full length. The design also supports a bacteria retentive air inlet which removes the possibility of air borne contamination and automatic shut off valve prevents air trapping in the fluid line once the measured volume chamber is empty and it also indicates the level of liquid.

The product also contains a 15µm fluid filter which reduces the risk of particulate matter.

The device is compatible with IV Cannula, Infusion Bottle and 3-way stopcock.

Product Description

Sterile, Single Use – Infusion Burette Set consist of transparent cylindrical calibrated measured volume chamber with printed graduation marking, Injection port for intermittent connection, Automatic Shut off valve which prevents air trapping in the fluid line once the measured volume chamber is empty and it also indicates the level of the liquid and 15µm fluid filter to reduce particulate matter.

Product Intended Use

Sterile, Single Use – Infusion Burette Set is used to administer Intravenous fluid and medicines into human circulating system by using Intravenous Catheter. The Burette Infusion Sets are designed to administer measured volume of infusion fluid through gravity method.

Product Features

- Clear and Transparent, Kink Resistance Tubing.
- Soft Transparent Cylindrical Shape with Calibrated measured volume chamber with sharp and bold marking scale.
- Injection port for intermittent medication.
- Automatic shutoff valve prevents air trapping in the fluid line once the measured volume chamber is empty and it also indicates the level of liquid.
- 15µm fluid filters reduces particulate matter.
- Efficient roller controller ensures better flow control.
- Luer locking system for secure connection.
- Free from Air and Liquid leakage
- Pyrogen Free
- Single Use, Disposable
- Sterile

Product Types and Variants

Burette Infusion Set

Size	Catalogue #
100ml	F33400010180
110ml	F33410011180
150ml	F33420015180



Our Company



Our passion is excellence, reliability, and creating innovative solutions.

Since our establishment in 1980, we have committed to ourselves to making the healthcare environment safe and clean by manufacturing the latest and highest-quality medical disposable.

We achieved that through a dedicated and integrated team that strives to make our company a leader in the field.

We have achieved many major accomplishments in the past, including owning the largest medical injection syringe factory in the Gulf and contributing to the fight against the coronavirus pandemic, by manufacturing over 60 million vaccine doses for Saudi Arabia and Dubai.

With many other promising projects and accomplishments that have made us a focus of attention and a hub for creative partners, we aim to be the first entity in our field, and we spare no effort in continuous development and innovation.

RAMY

SAUDI
MADE

